



O R M E S A[®]


design and quality made in Italy




M E A S U R E S




BIRILLO 1

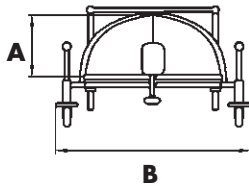

B

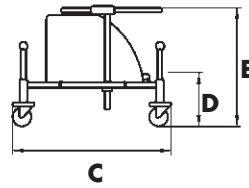

C

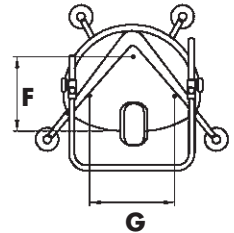

G

THE NUMBERS DIVIDED BY THE DASH SPECIFY A MINIMUM AND A MAXIMUM ADJUSTMENT

Sizes	Weight of the aid	Max. user mass	A	B	C	D	E	F	G
1	kg 8,5 Lbs. 18.7	kg 30 Lbs. 66	cm 18 7"	cm 59 23"	cm 50 19.6"	cm 19-27 7.5-9.8"	cm 36-50 14-19.7"	cm 18 7"	cm 22 8.6"
2	kg 10 Lbs. 22	kg 40 Lbs. 88	cm 19 7.4"	cm 65 25.6"	cm 51 20"	cm 22-30 8.6-11.8"	cm 40-56 15.7-22"	cm 20 7.8"	cm 26 10.2"
3	kg 11 Lbs. 24	kg 50 Lbs. 110	cm 21 8.3"	cm 73 28.7"	cm 59 23.2"	cm 22-34 8.6-13.4"	cm 42-62 16.5-24.4"	cm 23 9"	cm 30 11.8"
4	kg 12,8 Lbs. 28	kg 60 Lbs. 132	cm 28 11"	cm 78 30.7"	cm 60 23.6"	cm 22-40 8.6-15.7"	cm 50-76 19.6-30"	cm 30 11.8"	cm 38 15"

Reproduction, even partial, is forbidden. The manufacturer reserves the right to make any changes without prior notice.



BIRILLO 2



BIRILLO 3



BIRILLO 4

A C C E S S O R I E S



860 - concave headrest
adjustable in height



824 - transparent table

BIRILLO safeness is in compliance with the **93/42 European Directive**:

- The **frame**, adjustable in height, is in steel and it is painted with epoxy-atossic varnish;
- Fire-retardant **upholstery**;
- The **seat** is padded and covered with plastic material, resisting to wear and tear and cleaning;
- The elastic **wall-protector buffers** are over the wheels;
- The **solid tyres** are castored and they run on ball bearings;
- The **handlebar** is adjustable in height and depth.

BIRILLO was conceived and thought to let children from 2 to 8 years old precociously move by themselves in the sitting position. Its **upholstery's fire-retardant**. The **triangular sitting position** prevents from the extreme hip and trunk extension and lets the child have his head and trunk well lined up. According to this, the aid helps the child having a **right posture** and his **moving impulse**. The **abduction block** keeps the pelvis in the triangular sitting. Its **posture stimulates the movement** of the arms and the rotation of the trunk. Thanks to its **shape** and **colours**, **BIRILLO** is **accepted as a "toy"** by the child and this predisposes to a serene use, without fear, by **stimulating his curiosity of exploring the surrounding environment**. For this reason **BIRILLO** plays a very important part in the **development of the perceptive and motor visual functions of the child**, which are necessary to know and represent the space.



COMPANY
WITH QUALITY MANAGEMENT
SYSTEM CERTIFIED BY DNV
=ISO 9001:2000=



Authorized Dealer



S.r.l. - Via A. Da Sangallo, 1 - PACIANA
06034 FOLIGNO (PG) ITALY - TEL. +39 0742 22927 (r.a.) FAX +39 0742 22637
E-mail: info@ormesa.com - www.ormesa.com