











Freewheel or Fixed pinion

The **freewheel** allows the movement on the direction of travel only during the pedaling, so **it doesn't force the cyclist to undergo the movement on the legs**. With a freewheel-equipped bicycle, the legs can rest while the bike maintains its forward motion.

On the other hand, a **fixed pinion** setup causes the movement of the rear wheel to directly activate the pedal movement. As a result, it leads to **passive mobilization of the lower limbs**, even if they are not actively engaged in pedaling, e.g. while the caregiver pushes the bike.



Differential

The Differential, available as an option from the Small size, is a device that **greatly assists the user when navigating a curve**, allowing the two rear wheels to move at different speeds (the inner rear wheel covers a shorter distance compared to the outer one).

Upon request, the tricycle can be provided with either a **differential and freewheel** or a differential and fixed pinion.

The differential is included in the configuration with the three speed gear.



Three-speed gearbox

Three-speed gearbox gives the ability to very easily adjust effort and speed according to the tricyclist's ability.

With the gearbox you have the advantages of the freewheel version with the possibility of reverse gear, with coaster brake.

It can be activated while the triciclist is stationary. The display shows the current version of the gear.

Available from size Medium.

Sizes and colours

5 sizes and 4 colours to meet the needs of kids and a



dults













Anatomical supports

The supports offer posterior/lateral containment to the trunk and the pelvis and, thanks to the strap with closure, also good front containment.

Suitable for those users who need a help to keep the trunk in an erect and median position and stay on the bike more easily.

for the Trunk

Adjustable in height (via the 883 rod) in tilt and anterior-posterior direction via the joint.

for the Pelvis

Adjustable in height (via the 883 rod).









Front brake

Operated by a lever on the right or left grip handle. The distance between the brake lever and the grip handle is adjustable by screw.



Blocking system

Allows you to lock the tricycle, essential for getting on and off the tricycle. For sizes M, L and XL



Bell and ergonomic handles

For better grip and comfort.



Front Lamp

2 functions: fixed and flashing. For sizes M, I and XI



Back Lamp

With USB charging. 5 different light modes. With battery status indicator and usb cable. Road safety is granted also by reflective details along the basket



Mesh basket

Located on the back of the tricycle has an easily removable rigid base. The upholstery with reflective details is attached with snaps and is removable for washing.









Ergonomic saddle

Designed for maximum comfort. Adjustable in height, in tilt and anterior-posterior direction.




Mudguards

Rain, dust and mud are no problem. For sizes S, M, L and XL.



Handlebar

Adjustable in height, in tilt, anterior-posterior direction and reversible.

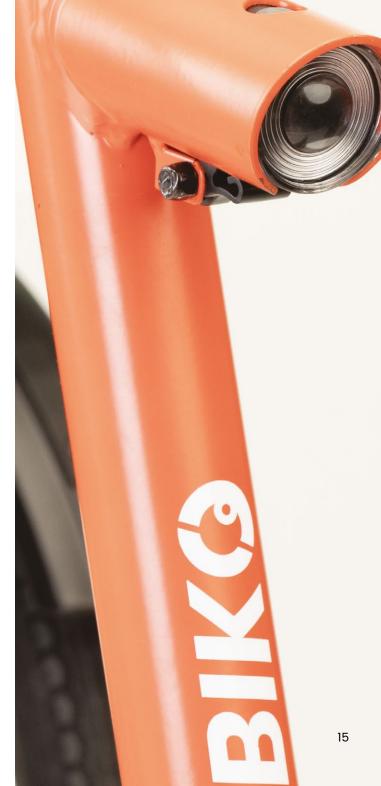





Components



840 Pedal cages with straps.To keep the forefoot positioned on the pedal.



841 Pedal cage with lenght and width adjustable foot imprints and fastening straps, with support kit.

It can be ordered per piece.



846 Ankle support adjustable in height with support kit.

(for use with the 841 pedal cage). To stabilize the tibiotarsal joint. It can be ordered per piece.



842P Pelvic support with belt*. Height adjustable thanks to 883



842T Multi-adjustable Trunk support with belt* and joint.

Adjustable in inclination, in anterior-posterior direction and in height thanks to 883



883 Rod adjustable in height for thoracic/pelvic support (single rod).

*Warning: It facilitates balance, but should not be used to support a user who doesn't have any trunk control

The postural components must be chosen and adjusted by a health professional who evaluates the residual abilities of the person; supervision of an adult, family member or qualified person, is always required during the use of the Tricycle.

Components



877 Flat backrest



843 Abductive saddle with holder for supports



adjustable glove independently.

Make a binding to the handlebar in case of an ineffective grip, it can also be used individually. With reflective details and easily removable for washing.



845 Forearm support It supports the forearm. It can be ordered per piece.



923 Additional assistant handleAllows the Care Giver to accompany the

child in riding the tricycle by assisting in pushing or slowing down, if necessary, without the ability to direct the handlebar. Easily removable for transportation.



878 MTB tires (3 pcs.) in alternative to standard tyres. Except size mini.



879 Front basketFor sizes mini, Small and Medium.





880 Adjustable ergonomic handles

They have the function and purpose of offering the user a grip closer to the body and more grip points than the standard handlebar. Useful for those users who have difficulty extending their elbows completely. They can be adjusted in height/inclination to accommodate various positions. They can be adjusted asymmetrically or independently, to accommodate any need for important asymmetry between the two upper limbs.



881 Steering limiter

It reduces the steering angle by preventing the steering from being inadvertently oriented in the transverse direction. For users who need help to manage the steering.



882 Steering limiter with spring

It reduces the steering angle by preventing the steering from being inadvertently oriented in the transverse direction. The return of the handlebar to the central position is facilitated by two springs that automatically orient the handlebar towards the forward-straight direction. For users who need help to manage the steering.

Components only choosable when ordering the tricycle



847 Cranks adjustable in length

For sizes M, L and XL. Decreases the distance between saddle and pedal of 3,5 cm and 5 cm/1,3 in and 1,9 in (see sizes pages).



864 Disc brake

applicable only on versions with differential and without gearbox.

Standard for sizes mini and Small.

sizes and weight

mini

45 cm~17.7 in

62 cm~24,4 in

39-53 cm

15,3-20,8 in

12 cm~4,7 in

103 cm~40,3 in

36 cm~14.2 in

90-110 cm

35,4-43,3 in

15,5 kg~34,1 lb

30 kg~66,1 lb

41x18x6,6 cm

16-7-2.6 in

15 kg~33 lb

A: handlebar width	
B: width encumbrance	
C1: pedal to saddle height	

D: ground to chain guard height

E: lenght encumbrance

F: seatpost to head tube distance

G: folded tricycle width

I: folded tricycle height

H: user height

Total weight

Maximum load

Storage basket LxMxN

Storage basket maximum load

* with standard saddle







Medium

55 cm~21.6 in

75 cm~29,5 in

59-72* cm

23,2-28,3 in

16 cm~6,3 in

136 cm~53,5 in

45 cm~17.7 in

81 cm~31.8 in

90 cm~35,4 in

120-145 cm

23 kg~50,7 lb

60 kg~132,2 lb

19.5-10.7-3.6 in

25 kg~55 lb

47,2-57 in

Small

50 cm~19,6 in

69 cm~27 in

47-63 cm

18,5-24,8 in

16 cm~6,3 in

41 cm~16 in

74 cm~29 in

80 cm~31.5 in

100-125 cm

39,3-49,2 in

20 kg~44 lb


45 kg~99,2 lb

17.6-9.3-2.6 in

25 kg~55 lb

44,9x23,6x6,6 cm

119 cm~46,8 in

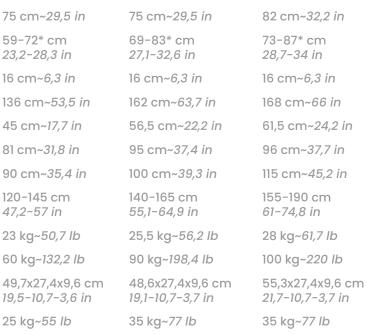


Large

60 cm~23,6 in

Storage basket view from above internal measures





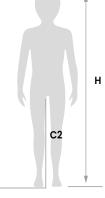

Extra Large

64 cm~25,2 in

C2= leg length







A special thanks to the families who participated with enthusiasm in our project. It is thanks to our customers and their stories that we find the energy and stimuli to evolve day by day.

(€ complies with MDR (EU) 2017/745

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 13485 =

Ormesa srl - Via delle Industrie 6/8 | 06034 Foligno (PG) - ITALY Tel. +39 0742 22927 | FAX +39 0742 22637 | info@ormesa.com

www.ormesa.com