

It stimulates and helps the step pattern. It trains the user to maintain the balance of the trunk and pelvis. It organizes a wider and structured vision of space, enhances the perception of distances and environment that surrounds the user. It develops physical self-sufficiency and psychological independence by experimenting the possibility of acting independently without the help of other people.

Grillo compensates the user lack of balance thanks to its ergonomic support and wide base, it allows the user to concentrate his/her residual abilities on walking.

The wide range of components and supports offers maximum versatility. All accessories are interchangeable quickly and without the use of tools to suit different clinical pictures.



Product registered FDA (USA)

Product marked CE


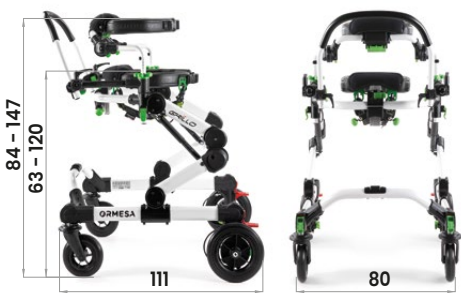


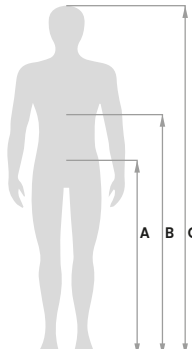


- **LIGHT ALUMINUM FRAME**, foldable for an easy transportation with PATENTED HEIGHT ADJUSTMENT SYSTEM.
- **PELVIC SUPPORT and THORACIC SUPPORT**, padded, adjustable in height, width and circumference. They allow an easy access. They hold the user just like two HANDS giving him/her a sense of security and stability allowing the user to concentrate on walking without any fear. The Thoracic support is adjustable in tilt and backward-forward.
- **ERGONOMIC HARNESS**, comfortable, adjustable and easily removable for washing. It supports the user in case of necessity without generating painful pressures.
- **HANDLEBAR**, removable and adjustable in height, it is used as a reference and as an armrest.

- **ARM SUPPORTS**, adjustable in height, tilting, rotation and depth.
- Exclusive **ORMESA DESIGN WHEELS**, light and puncture proof : 2 front swivelling, Ø15 cm (sizes mini and Small) and Ø 17.5 cm (sizes Medium and Large) and 2 rear fixed, Ø 20 cm (sizes mini and Small) and Ø 25 cm (sizes Medium and Large).
- **DIRECTIONAL LOCK** on front wheels, **ADJUSTABLE DRAG**, **ANTI-ROLLBACK DEVICE** and **WHEELS STOP** on rear wheels, easy to operate by the red pedal
- **SEMIELASTIC BUMPERS**.



Attention!
Read carefully the user and maintenance manual

FRONT DRIVE	REAR DRIVE	USER DIMENSIONS:
<p>Inner circumference pelvic and thoracic support</p>  <p>from cm 81 to cm 140</p> 	<p>Inner circumference pelvic and thoracic support</p>  <p>from cm 81 to cm 140</p> 	<p>USER HEIGHT (C): cm 150-190</p> 
<p>WEIGHT OF THE AID kg. 29,5 MAX LOAD kg. 110</p>	<p>WEIGHT OF THE AID kg. 28,5 MAX LOAD kg. 110</p>	<p>IMPORTANT: Make sure that measures A and B indicated in the figure do not exceed the suggested range indicated for the supports.</p>

download and send the form to
info@ormesa.com

DISTRIBUTOR/AGENT:	DATE OF ORDER:
	DATE OF DELIVERY:
	NAME OF THE ADDRESSEE/USER:

ORDER FORM

GRILLO

SIZE **L**

Front Drive

Rear Drive

CHOOSE THE CONFIGURATION >	Front Drive			Rear Drive		
	PT	PA	P	PT	PA	P
ADDITIONAL COMPONENTS ✓ 925, 924, 930: you can choose other size	<input type="checkbox"/> Including: > thoracic support > pelvic support > ergonomic harness > handlebar	<input type="checkbox"/> Including: > pelvic support > ergonomic harness > arm supports	<input type="checkbox"/> Including: > pelvic support > ergonomic harness > handlebar	<input type="checkbox"/> Including: > thoracic support > pelvic support > ergonomic harness > handles	<input type="checkbox"/> Including: > pelvic support > ergonomic harness > arm supports	<input type="checkbox"/> Including: > pelvic support > ergonomic harness > handles
865 Multiadjustable headrest	NA	NA	NA	<input type="checkbox"/> available in combination with thoracic support 925	<input type="checkbox"/> available in combination with thoracic support 925	<input type="checkbox"/> available in combination with thoracic support 925
925 Thoracic support	●	<input type="checkbox"/>	<input type="checkbox"/>	●	<input type="checkbox"/>	<input type="checkbox"/>
924 Pelvic support	●	●	●	●	●	●
930 Ergonomic harness	●	●	●	●	●	●
926 Adjustable handlebar	●	<input type="checkbox"/>	●	NA	NA	NA
809 Arm supports	<input type="checkbox"/>	●	<input type="checkbox"/>	<input type="checkbox"/> with joints	●	<input type="checkbox"/> with joints
943 Arm straps	<input type="checkbox"/> available in combination with arm supports 809	<input type="checkbox"/>	<input type="checkbox"/> available in combination with arm supports 809	<input type="checkbox"/> available in combination with arm supports 809	<input type="checkbox"/>	<input type="checkbox"/> available in combination with arm supports 809
939 Ergonomic handles	<input type="checkbox"/>	●	<input type="checkbox"/>	<input type="checkbox"/> with joints	●	<input type="checkbox"/> with joints
927 Removable knobs	NA	NA	NA	●	<input type="checkbox"/>	●
815 Proximal abductor	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
946 Ergonomic saddle	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
890-C Leg divider with thigh loops	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
890-SC Distal abductor	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
810 Weighting bars	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
928 Rear narrow base wheels with parking brake pedals. ONLY FOR INDOOR USE	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>