

The alternative to being bedridden



Marked

Tie down hooks (4 rings) for the front facing transport of the user in moving vehicles (private cars, conformed buses etc.)
AS STANDARD



Attention!
Read carefully the user and maintenance manual

DISTRIBUTOR/AGENT:	DATE OF ORDER:
	DATE OF DELIVERY:
	NAME OF THE ADDRESSEE/USER:

FEATURES

- **FRAME** made of chromed steel, painted with non-toxic epoxy powder paint. Narrow base to facilitate entering elevators.
- **The seat and backrest are covered with visco-elastic fire-proof material. It is resistant, breathable and anti-slip. It is made of a special and exclusive ORMESA formula:**
 - 1) It adapts to the user's body shapes and it prevents its sliding, with the result of an **EXCELLENT COMFORT and KEEPING POSTURE**, even when changing the backrest recline and tilt (it adapts to the user's body shapes maintaining its support function).
 - 2) Easily washable to keep the **HYGIENE** of the user; it lets **AIR PASS** keeping the **OPTIMAL MICROCLIMATE** that **HELPS PREVENTING BED SORES**.
 - 3) Durable and fireproof, **GUARANTEE OF DURABILITY and SAFETY**.
 - 4) Ergonomic, breathable, non-toxic, fireproof **UPHOLSTERY**; It does not interfere with the benefits of the visco-elastic cover. The upholstery is fixed with buttons, easy to remove for machine washing.
- Ergonomic **BACKREST**, it holds the shapes of the user. **The special structure of the backrest gives a proper alignment to the trunk and it doesn't generate any pressure on the spine.** Adjustable in tilt by levers with safety pin. The pin allows blocking the levers to maintain the proper adjustment for the user, height adjustable, folding for transportation.
- Low **SEAT** and slightly inclined. Adjustable in depth on four positions.
- Wide and wrap **HEADREST**, padded, multi adjustable (in height, tilt, depth, lateral traslation); thus assuring the **proper alignment between head and trunk**, without interfering with the user's view. Anti-bacterial and Hygienic.
- **ARMRESTS**, padded, easily and stably adjustable in height thanks to the lever.
- **PUSH HANDLE** adjustable in tilt/ height, it helps the care.
- **CONTINUOUS AND WITHOUT EFFORT ADJUSTMENTS**: tilting, backrest recline and legrest tilt are operated by independent gas springs.
- **SIDES** with scissor opening. Thanks to the special joint, the armrests are positioned behind the backrest, eliminating completely all obstacles to the transfer.
- **LEG RESTS**: ergonomic, removable, turnable and adjustable in tilt without effort, activated with gas springs.
- **CALFRESTS** adjustable in depth and free tilt to accommodate the movements of the leg.
- **FOOT RESTS**, adjustable in flexion-extension and height.
- **8" front solid wheels and 22" rear wheels with brake.** Anti-tip device and pedal to climb steps.


Please check off the interested products:

COLORS MEASURES

FRAME COLOUR: white black

UPHOLSTERY COLOUR: dark grey light grey

SIZES: 36 40 45 50

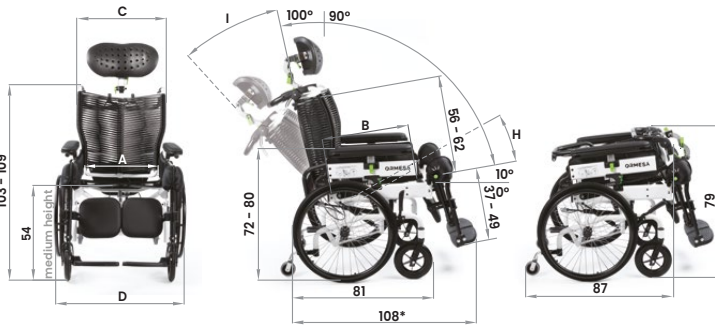

*footrest adjusted to the maximum length

Juditta B60	size 36	size 40	size 45	size 50
A seat bottom width	cm 36	cm 40	cm 45	cm 50
B seat bottom depth	cm 40-43-46-49	cm 40-43-46-49	cm 40-43-46-49	cm 40-43-46-49
C backrest width	cm 47	cm 47	cm 52	cm 57
D overall width	cm 69	cm 69	cm 74	cm 79
H seat bottom inclination-tilt	from 10° to 35°	from 10° to 35°	from 10° to 35°	from 10° to 35°
I backrest inclination	from 10° to 38°	from 10° to 38°	from 10° to 38°	from 10° to 38°
Total weight	kg 40	kg 40	kg 42	kg 43
Max capacity	kg 120	kg 120	kg 130	kg 140

Hyphen separated numbers indicate a minimum and maximum range of adjustment

ADDITIONAL COMPONENTS

824-R Reclinable tray, adjustable in tilt, backwards-forwards, with security block

It is made by resistant and hygienic material, design with no edges, it doesn't keep residual food.



894 45° Pelvic belt



834 Padded abduction block



892 Footplate covering



905 Set of hand brakes with single lever



914 Footrest lock




896 Tray 824 covering



916 Padded leg panel



936 Headrest cover

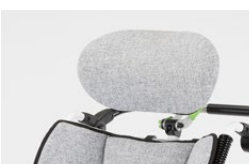

938 Water-proof cover

