

DISTRIBUTOR/AGENT:	DATE OF ORDER:
	DATE OF DELIVERY:
	NAME OF THE ADDRESSEE/USER:

FEATURES

- **FRAME** made of chromed steel, painted with non-toxic epoxy powder paint. Narrow base to facilitate entering elevators.
- **The seat and backrest are covered with visco-elastic fire-proof material. It is resistant, breathable and anti-slip. It is made of a special and exclusive ORMESA formula:**
 - 1) It adapts to the user's body shapes and it prevents its sliding, with the result of an **EXCELLENT COMFORT and KEEPING POSTURE**, even when changing the backrest recline and tilt (it adapts to the user's body shapes maintaining its support function).
 - 2) Easily washable to keep the **HYGIENE** of the user; it lets **AIR PASS** keeping the **OPTIMAL MICROCLIMATE** that **HELPS PREVENTING BED SORES**.
 - 3) Durable and fireproof, **GUARANTEE OF DURABILITY and SAFETY**.
- Ergonomic, breathable, non-toxic, fireproof **UPHOLSTERY**; It does not interfere with the benefits of the visco-elastic cover. The upholstery is fixed with buttons, easy to remove for machine washing.
- Ergonomic **BACKREST**, it holds the shapes of the user. **The special structure of the backrest gives a proper alignment to the trunk and it doesn't generate any pressure on the spine.** Adjustable in tilt by levers with safety pin. The pin allows blocking the levers to maintain the proper adjustment for the user, height adjustable, folding for transportation.
- Low **SEAT** and slightly inclined. Adjustable in depth on four positions.
- Wide and wrap **HEADREST**, padded, multi adjustable (in height, tilt, depth, lateral traslation); thus assuring the **proper alignment between head and trunk**, without interfering with the user's view. Anti-bacterial and Hygienic.
- **ARMRESTS**, padded, easily and stably adjustable in height thanks to the lever.
- **PUSH HANDLE** adjustable in tilt/ height, it helps the care.
- **CONTINUOUS AND WITHOUT EFFORT ADJUSTMENTS**: tilting, backrest recline and legrest tilt are operated by independent gas springs.
- **SIDES** with scissor opening. Thanks to the special joint, the armrests are positioned behind the backrest, eliminating completely all obstacles to the transfer.
- **LEG RESTS**: ergonomic, removable, turnable and adjustable in tilt without effort, activated with gas springs.
- **CALFRETS** adjustable in depth and free tilt to accommodate the movements of the leg.
- **FOOT RESTS**, adjustable in flexion-extension and height.
- **8" front solid swiveling wheels and 12" rear wheels**; drum brakes are operated by a single control pedal.



Marked **CE**

Tie down hooks (4 rings) for the front facing transport of the user in moving vehicles (private cars, conformed buses etc.)

AS STANDARD



Please check off the interested products:



Attention!
Read carefully the user and maintenance manual

SIZE




UPHOLSTERY COLOUR:



FRAME COLOUR:



MEASURES

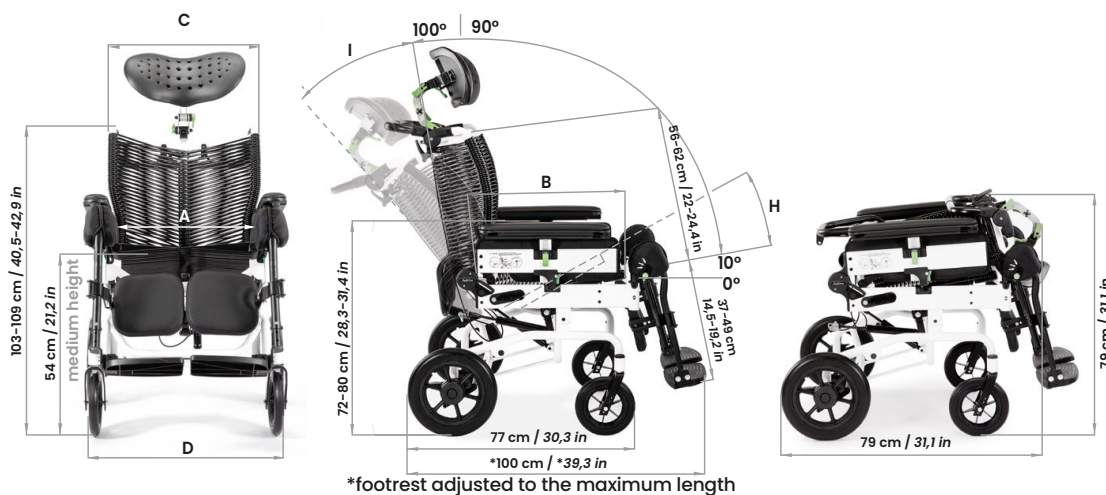

*footrest adjusted to the maximum length

	A seat bottom width	B seat bottom depth	C backrest width	D overall width	H seat bottom inclination-tilt	I backrest inclination	Total weight	Max. Capacity	User height
size 36	36 cm 14,1 in	40-43-46-49 cm 15,7-16,9-18,1-19,2 in	47 cm 18,5 in	62 cm 24,4 in	from 10° to 35°	from 10° to 38°	38 kg 83,7 lb	120 kg 264,5 lb	from 140 to 190 cm from 55,1 to 74,8 in
size 40	40 cm 15,7 in	40-43-46-49 cm 15,7-16,9-18,1-19,2 in	47 cm 18,5 in	62 cm 24,4 in	from 10° to 35°	from 10° to 38°	38 kg 83,7 lb	120 kg 264,5 lb	from 140 to 190 cm from 55,1 to 74,8 in
size 45	45 cm 17,7 in	40-43-46-49 cm 15,7-16,9-18,1-19,2 in	52 cm 20,4 in	67 cm 26,3 in	from 10° to 35°	from 10° to 38°	41 kg 90,3 lb	130 kg 286,6 lb	from 140 to 190 cm from 55,1 to 74,8 in
size 50	50 cm 19,6 in	40-43-46-49 cm 15,7-16,9-18,1-19,2 in	57 cm 22,4 in	72 cm 28,3 in	from 10° to 35°	from 10° to 38°	42 kg 92,5 lb	140 kg 308,6 lb	from 140 to 190 cm from 55,1 to 74,8 in

Hyphen separated numbers indicate a minimum and maximum range of adjustment

ADDITIONAL COMPONENTS

- ☐ **824-R**
Reclinable tray,
adjustable in tilt,
backwards -
forwards, with
security block

It is made by resistant
and hygienic material,
design with no edges,
it doesn't keep
residual food.



- ☐ **894** 45° Pelvic belt




- ☐ **834** Padded
abduction block



- ☐ **892** Footplate
covering



- ☐ **905** Seat of hand
brakes with single lever

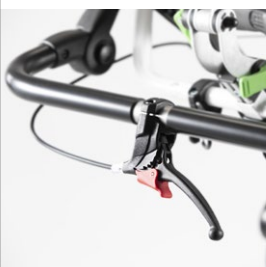

- ☐ **914** Footrest lock



- ☐ **896** Tray covering



- ☐ **916** Padded leg
panel



- ☐ **932** Anti-tip device



- ☐ **933** IV Pole



- ☐ **936** Headrest
cover



- ☐ **938**
Water-proof cover



- 975** Multi adjustable
side-temporo support. It
can be positioned on the right or left side.

- ☐ 1 pc. ☐ 2 pc.



- ☐ **976** Multi adjustable
side-temporo-fronto
support. It can be positioned on
the right or left side.



Notes: _____
