

INGLESE - cod.220615 - 22.04.26



**ORMESA®**  
• MADE IN ITALY SINCE 1980 •

**GRILLO**  
*adaptive seating*



# GRILLO

*adaptive seating*

## More than just a high chair.

**Grillo Adaptive Seating** by Ormesa is a **compact activity chair** with adjustable configuration for children and teens, ideal to be used in many contexts, such as **carrying out daily activities** at school and home.

The **international height adjustment patented system** typical of the “**GRILLO product family**” allows for a **compact, space-saving base** while offering a wide range of height adjustments.

**Grillo Adaptive Seating** has a **extremely breathable, durable, easily washable** and **fire retardant** load-bearing upholstery.

The removable cover is made with **highly innovative padding**: the material allows an **optimal weight distribution**, very **high breathability**, and has an **excellent comfort**. Available in **three colors: red, pink, and blue**.

With a **wide range of available components**, it easily adapts to any need. Available in **three sizes**.

PRODUCT  
MARKED 

GRILLO ADAPTIVE SEATING IS A MEDICAL DEVICE COMPLYING WITH THE SAFETY REQUIREMENTS DEFINED BY REGULATION (EU) 2017/745



# Three configurations to meet every need

**Grillo Adaptive Seating** by Ormesa is an assistive device **designed to meet the needs of care givers and children**: parents know that the activity chair becomes a fundamental tool in the development of the child, for the performance of any activity; the child implements the ability to see the surrounding world with a new perspective, to feel more involved in everyday activities and to communicate with others in a whole new way. Three different configurations are created to meet every need:

- › Grillo Adaptive Seating with **component 863 external postural headrest** suitable for children who need an adjustable support in inclination and depth.
- › Grillo Adaptive Seating with **component 865 padded headrest** or **component 835 shaped headrest**, depending on the need for lateral and occipital head restraint.





**863 Postural headrest**

That adapts to the conformation of the head. Adjustable in height, tilt, and depth.



**865 Padded headrest**  
adjustable in height



**835 Shaped headrest**  
with lateral and occipital support, adjustable in height.

# Standard and lower base



## The versatility you need

Grillo Adaptive Seating has an **aluminum folding frame**, painted with non-toxic epoxy powders.

The **patented gas-spring-operated height adjustment system** allows to change the height of the chair with the child seated: ideal for getting him or her to access to different table tops and work surfaces with ease and comfort, and to bring him or her closer to the floor for easy descent.

**The safety system avoids the inadvertent activation.**

Thanks to its **smart design, it adapts to its surroundings**, allowing the child to participate in daily activities and interact with family and other children.



**THE 4 SWIVELLING WHEELS, with the NEW rubber coating** designed to prevent dirt from getting trapped in them and all equipped with parking brake.

They allow you to move the Grillo Adaptive Seating effortlessly while the child is sitting.




### Lower base

to get closer to the ground than is possible with the standard indoor base (*useful in community settings such as nursery and primary schools.*) It keeps the same features as the standard base and it is available for sizes mini and small.



# High performance

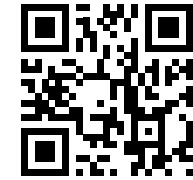

## Cutting-edge double cover



The net and non-removable supporting structure, in combination with the technological upholstery, besides offering **excellent torso and pelvis stability**, guarantees a **very high level of breathability**, to the benefit of comfort and health of the skin.



- **Optimal "self-supporting effect"**: the highly innovative padding structure allows **optimal weight distribution** and has an excellent comfort. Pressure is optimally distributed with a significant reduction in point pressure.
- It has a **compressive stress resistance** (Compression stress value) (DIN EN ISO 3386) higher than traditional foam rubber.
- **Lightweight and breathable**: provides extremely high air circulation, **reducing sweating** and promoting thermoregulation.
- **Hygienic**: easy to wash with fast drying times.

- **Abrasion-resistant**: The high-strength coating texture maintains the performance even under high loads (ISO 12947 test).
- **Fire-retardant**: According to EN 1021-1:2014 and EN 1021-2:2014..
- **Oeko-Tex certificate**: for all fabrics that are in direct contact with the skin, as a guarantee of quality and health safety.

# Technical features and adjustments



## Backrest

Reclines with gas spring easily and securely with physiological fulcrum, to make the adjustment comfortable. Adjustment in height by sliding, without any risk of losing the screw. Foldable for transportation.



## Seat

Ergonomic and tilting with gas spring for greater comfort. Adjustable in depth by sliding, without any risk of losing the screw, and in width by multi-adjustable padded pelvic side supports.



## Multi-adjustable pelvic side supports

For greater containment. Adjustable in rotation, width and depth.



**NEW**

## Multi-adjustable armrests

Soft and adjustable in height, inclination, rotation and length.

## Legrest

Adjustable in height and tilt, in depth along with the seat with physiological fulcrum, to make the adjustment comfortable.




## Calf rests

Soft and adjustable in height and depth.

## Footrests

Tilting and adjustable in flexion-extension.

# Compact and stable adjustments without barriers

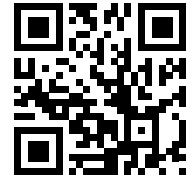

# Components



## **947 Pelvic belt with variable angle**

Available in three sizes. Provided as standard.



## **920 Four point pelvic belt**

Available in two sizes.



## **838 Trunk side supports**

Adjustable in height, width and rotation. Available in three sizes.

*Covered by breathable padding fabric, removable for washing.*



## **868 Wrappable and flexible trunk supports**

Adjustable in height, width and rotation. Available in two sizes.

*Covered by breathable padding fabric, removable for washing.*



## **961 Elbow side supports**

*Covered by breathable padding fabric, removable for washing.*



## **853 Vest harness**

Available in two sizes.

# Components



**853 slim Four point shaped harness**



**834R Adjustable abduction block**  
in depth and abduction. Available in two sizes.



**834N Narrow abduction block**  
Adjustable in depth along with the seat. Available in two sizes.



**824 Tray with wrap-around recess**  
Adjustable in depth, height and tilt with armrests.

*It is made of resistant and hygiene material, no edge design, it does not retain any residue of food.*




**960 Heel rests**



**962 Footrests + Heelrests**

# Grillo A.Seating

## with standard base: measurements (cm/In)

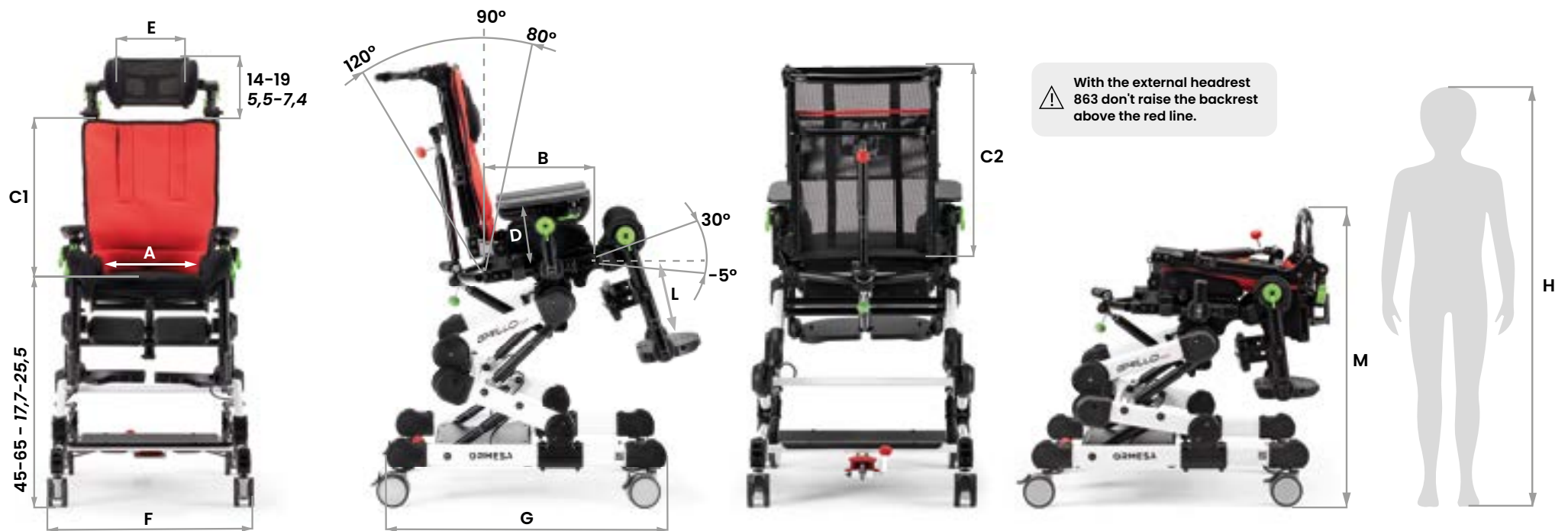


	A	B	C1	C2	D	E	F	G	L	M	user height (H)	total weight	max. load
<b>mini</b>	16-26 6,2-10,2	16-26 6,2-10,2	30-41 11,8-16,1	30-50 11,8-19,6	12-15 4,7-5,9	15 5,9	49 19,2	59 23,2	15-25 5,9-9,8	68 26,7	75-105 29,5-41,3	20,4kg 44,9 lb	25 kg 55,1 lb
<b>Small</b>	20-30 7,8-11,8	20-30 7,8-11,8	35-47 13,7-18,5	35-60 13,7-23,6	15-18 5,9-7	17 6,6	53 20,8	65 25,5	20-31 7,8-12,2	68 26,7	90-120 35,4-47,2	22,6 kg 49,8 lb	35 kg 77,1 lb
<b>Medium</b>	25-35 9,8-13,7	28-43 11-16,9	40-55 15,7-21,6	40-70 15,7-27,5	16-20 6,2-7,8	20 7,8	58 22,8	74 29,1	25-40 9,8-15,7	70 27,5	110-145 43,3-57	24,3 kg 53,5 lb	45 kg 99,2 lb

Hypen-separated numbers indicate the minimum and maximum range of adjustment

# Grillo A.Seating

## with lower base: measurements (cm/In)

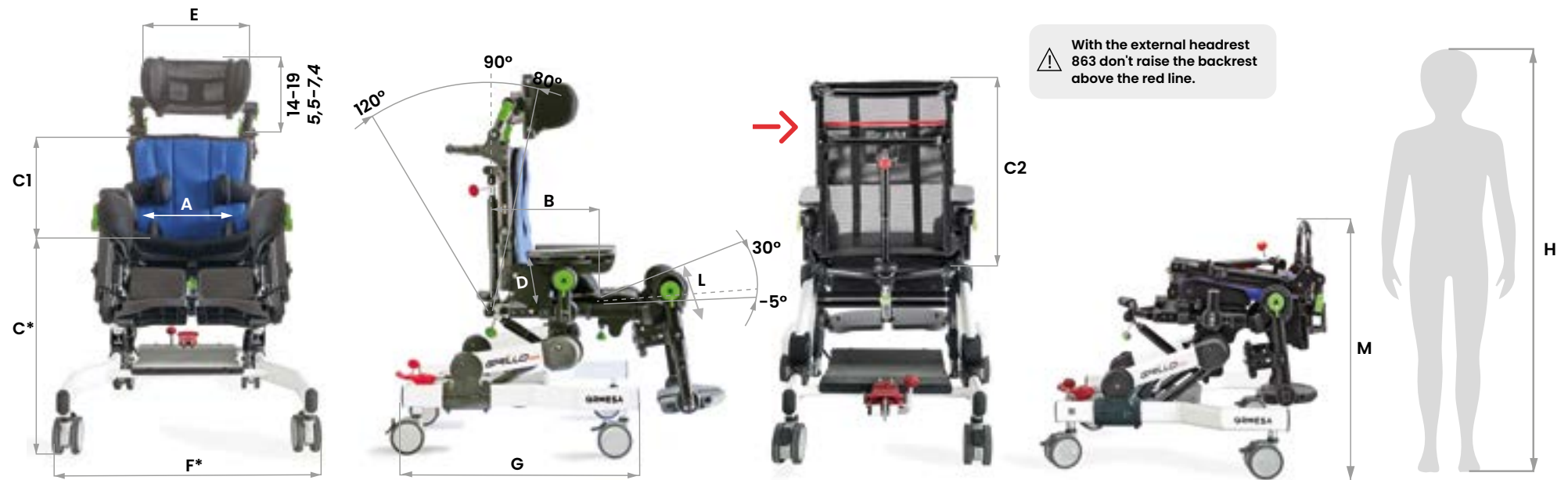


	A	B	C*	C1	C2	D	E	F*	G	L	M	User height (H)	Tot. weight	Max load
mini	16-26 6,2-10,2	16-26 6,2-10,2	30-40 11,8-15,7	30-41 11,8-16,1	30-50 11,8-19,6	12-15 4,7-5,9	15 5,9	55,5 21,8	57 22,4	15-25 5,9-9,8	58 22,8	75-105 29,5-41,3	18 kg 39,6 lb	25 kg 55,1 lb
Small	20-30 7,8-11,8	20-30 7,8-11,8	36-48 14-18,9	35-47 13,7-18,5	35-60 13,7-23,6	15-18 5,9-7	17 6,6	60 23,6	63 24,8	20-31 7,8-12,2	60 23,6	90-120 35,4-47,2	18,5 kg 40,7 lb	35 kg 77,1 lb

Hyphen-separated numbers indicate the minimum and maximum range of adjustment

Tessa and her story await you in the Ormesa official website, together with other protagonists of our fantastic global community.



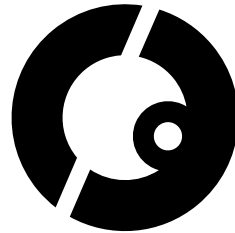

**A special thanks to the families who participated with enthusiasm in our project. It is thanks to our customers and their stories that we find the energy and stimuli to evolve day by day.**

**CE** complies with MDR (EU) 2017/745

COMPANY WITH  
QUALITY SYSTEM  
CERTIFIED BY DNV  
ISO 13485

---

**Ormesa srl** - Via delle Industrie 6/8 | 06034 Foligno (PG) - ITALY  
Tel. +39 0742 22927 | FAX +39 0742 22637 | info@ormesa.com

**www.ormesa.com**

